

[54] MASS SPECTROMETER SYSTEM AND METHOD FOR MATRIX-ASSISTED LASER DESORPTION MEASUREMENTS

[76] Inventor: Marvin L. Vestal, 7851 Chinon Cir., Houston, Tex. 77071

[21] Appl. No.: 278,405

[22] Filed: Jul. 21, 1994

[51] Int. Cl.⁶ G01N 24/00; G01N 35/02; H01J 49/04

[52] U.S. Cl. 436/47; 250/288; 250/289; 250/423 P; 250/442.11; 436/173; 436/181

[58] Field of Search 250/288, 289, 250/423 P, 442.11; 436/47, 173, 181

[56] References Cited

U.S. PATENT DOCUMENTS

4,076,982	2/1978	Ritter et al.	250/288
4,296,322	10/1981	Wechsung	250/423 P X
4,330,208	5/1982	Eloy	250/288 X
4,405,860	9/1983	Brunnee et al.	250/289 X
4,620,103	10/1986	Kambara et al.	250/288 X
4,740,692	4/1988	Yamamoto et al.	250/423 P X
4,988,879	1/1991	Zave et al.	250/423 P
5,041,725	8/1991	Izumi	250/288 X
5,045,694	9/1991	Beavis et al.	250/287
5,083,020	1/1992	Takahashi	250/288
5,160,840	11/1992	Vestal	250/282 X
5,288,644	2/1994	Beavis et al.	250/282 X
5,382,793	1/1995	Weinberger et al.	250/288

OTHER PUBLICATIONS

Article, "Metastable Decay of Peptides and Proteins in Matrix-Assisted Laser-Desorption Mass Spectrometry", Rapid Communications in Mass Spectrometry, vol. 5, 198-202 (1991).

Primary Examiner—Jeffrey R. Snay


Assistant Examiner—Arlen Soderquist

Attorney, Agent, or Firm—Browning, Bushman, Anderson & Brookhart

ABSTRACT

[57] The system for analyzing multiple samples includes a plurality of portable of sample supports each for accommodating a plurality of samples thereon, and an identification mechanism for identifying each sample location on each of the plurality of sample supports. The mass spectrometer is provided for analyzing each of the plurality of samples when positioned within a sample receiving chamber, and a laser source strikes each sample with a laser pulse to desorb and ionize sample molecules. The support transport mechanism provided for automatically inputting and outputting each of the sample supports from the sample receiving chamber of the mass spectrometer. A vacuum lock chamber receives the sample supports and maintains at least one of the sample supports within a controlled environment while samples on another of the plurality of sample supports are being struck with laser pulses. The computer is provided for recording test data from the mass spectrometer and for controlling the operation of the system.

74 Claims, 5 Drawing Sheets

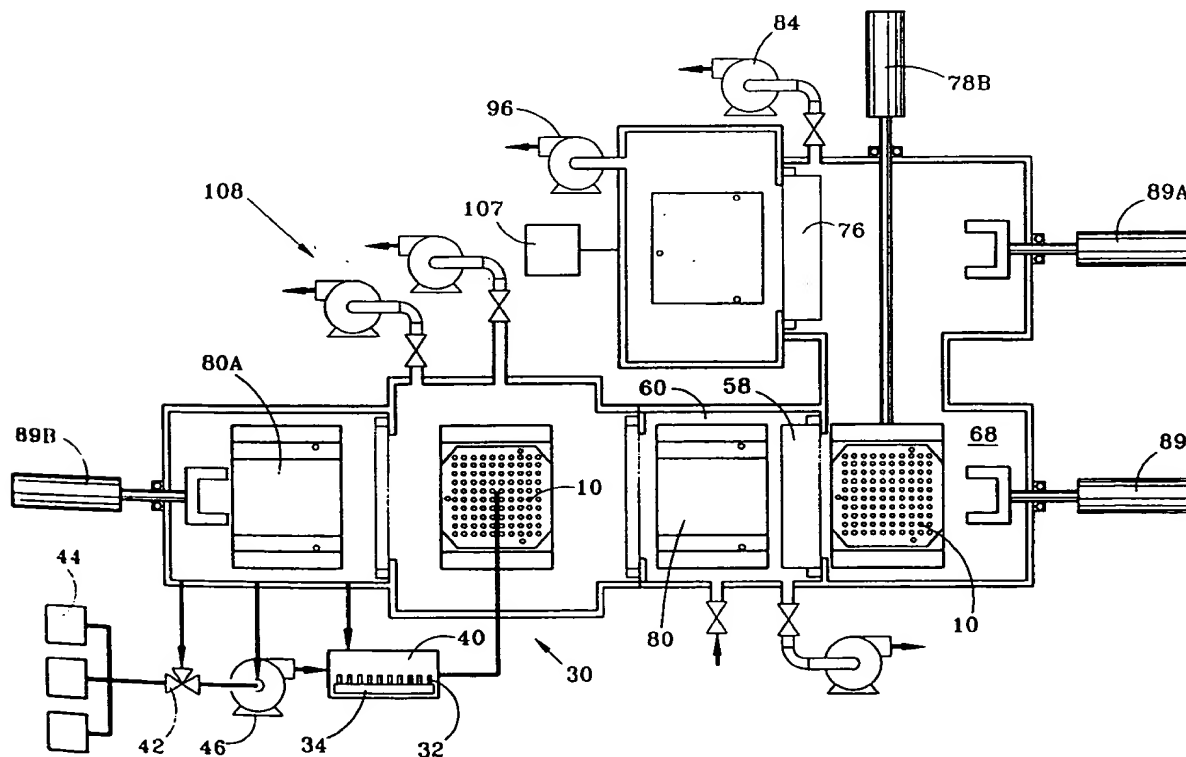

tenance fees
ears and six
reafter upon
ming of the
ment of the
or before the
l expire as of